**Jordan\***

CLEAN SMILE

Electric Toothbrush

Instruction manual  
Bruksanvisning  
Brugsanvisning  
Instruktionsbok  
Käyttöohje  
Instrukcja obsługi  
Handleiding  
Bedienungsanleitung  
Manuel d'instructions

TB-100GY

## **JORDAN CLEAN SMILE ELECTRIC TOOTHBRUSH**

The Jordan Clean Smile series is a new generation of electric toothbrushes that through SoftClean™ Technology delivers the right combination of oral hygiene features. Inspired by a Scandinavian design approach for simplicity and user-friendliness, the series combines the effectiveness of an electric toothbrush with the ease-of-use of a manual one.

---

### **SAFETY**

Please note the following instructions when using the appliance:

- The toothbrush can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should always be supervised when cleaning, and should not play with the appliance.

### **IMPORTANT INFO:**

- If the cord is damaged, it must be replaced and not used.
- If you are undergoing treatment for oral care, consult your dental professional prior to use.
- The toothbrush can only be charged by the charging unit provided with the appliance.
- When charging please make sure the toothbrush or charging unit does not come in contact with water or any other liquid.

### **OPERATION TIME**

Maximum operation time is approximately 80 minutes / 20 days when fully charged. The battery indicator on the handle will start blinking from 20 % battery capacity. Please recharge as soon as possible to preserve the life time of the battery. If however you continue to brush you might experience that the toothbrush might stop automatically.

### **CHARGING THE UNIT**

The toothbrush can be charged with the supplied USB cable.

**IMPORTANT:** Remember to close the lid in the bottom of the toothbrush after charging.

While charging the indicator will blink and will light up WHITE when the batteries are fully charged (2-3 hours). You can also connect the USB cable to your computer directly for charging. The power unit may warm up somewhat during recharging. This has no significance on the unit's operation. Over-charging is impossible.


### **BEFORE USE**

Before using your Jordan Clean Smile

toothbrush for the first time, charge the unit fully. After the power unit has been fully charged, attach the brush head to the power unit. The brush head must face the same direction as the on/off button. The toothbrush will stop if you turn it on while it is charging.

### **PRESSURE SENSOR**

When excessive pressure is applied when brushing, the pressure sensor indicator will appear on LED display and the oscillating speed will be reduced automatically to remind you to lower your pressure and also protecting your teeth and gum.
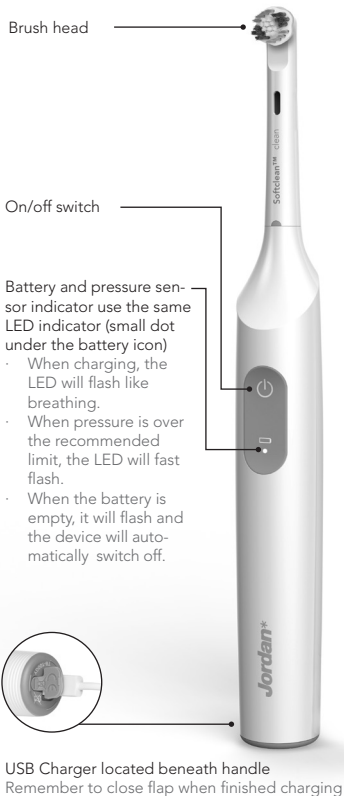
## **BATTERY**

This toothbrush uses a Lithium rechargeable battery. At the end of the toothbrush life time, you need to deliver the toothbrush for proper disposal according to your local authority instruction for environmental protection disposal.

## **PROTECTING THE ENVIRONMENT. CORRECT DISPOSAL OF THIS PRODUCT.**

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

## **POWER SUPPLY**

### **Jordan USB cable**

#### **Appliance:**

Lithium-ion battery 3.7V, charged by DC5V, 1A.

Distributed by Wilfa:

#### **WILFA AS**

Industriveien 25,  
1481 Hagan  
Norway  
[www.wilfa.com](http://www.wilfa.com)

#### **SUPPORT:**

[www.wilfa.com/support](http://www.wilfa.com/support)

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